# Dose Distribution of Gamma Irradiators

Seungwoo Park, Sang Hun Shin, Kihong Son, Chang Yeol Lee, Kum Bae Kim, Haijo Jung, Young Hoon Ji

Korea Institute of Radiological and Medical Sciences, Seoul, Republic of Korea E-mail: swpark@kirams.re.kr

keyword : Dose distribution, Gammacell 3000 Elan, BioBeam 8000

## Introduction

Gamma irradiator using Cs-137 have been widely utilized to the irradiation of cell, blood, and animal, and the dose measurement and education. The Gammacell 3000 Elan (Nordion International, Kanata, Ontario, Canada) irradiator was installed in 2003 with Cs-137 and dose rate of 3.2 Gy/min. And the BioBeam 8000 (Gamma-Service Medical GmbH, Leipzig, Germany) irradiator was installed in 2008 with Cs-137 and dose rate of 3.5 Gy/min.

Our purpose was to evaluate the practical associated dosimetric problems with inhomogeneous dose distribution within the irradiated volume in open air state using glass dosimeter and Gafchromic EBT film dosimeter[1-2] for routine Gamma irradiator dosimetry applications at the KIRAMS and the measurements were compared with each other. In addition, an user guideline for useful utilization of the device based on practical dosimetry will be prepared.

#### **Materials and Methods**

Dose distribution measurements were achieved by means of the glass rod dosimeter, and the Gafchromic EBT film. Irradiations were performed to a nominal absorbed dose of 0.5<sup>6</sup> Gy in commercially available two gamma irradiators (Gammacell 3000 Elan, and BioBeam 8000).

The uniformity dose distributions inside each irradiation chamber were determined by using the Gafchromic EBT film. The dosimetry of the Gafchromic EBT film was performed with irradiation dose not exceeding 7Gv in consideration of the characteristic of the Gafchromic EBT film. The glass dosimeters were located in the center and edge of the container and the Gafchromic EBT film were located in the center of the container. Reader used to reading of the glass dosimeter was to FGD-1000 Ace(Asahi Techno Glass Corporation, Shizuoka, Japan). Irradiated films were scanned and analyzed using EPSON Pro 1680 Expression scanner (Seiko Epson. Corporation, Nagano, Japan) in red channel.

#### Results

Minimum and maximum dose rate in the Gammacell 3000 irradiator of the glass dosimeters at central axis was measured to 4.06 Gy/min at the highest position, and 5.33 Gy/min at the middle position, respectively. Minimum and maximum dose rate in the BioBeam 8000 irradiator of the glass dosimeters at central axis was measured to 2.99 Gy/min at the highest position, and 3.66 Gy/min at the middle position,

respectively.

Fig. 1, Fig. 2 shows percent distributions based on reference value with the minimum dose rate position in order to compare dose uniformity within the irradiator. Maximum variation at each irradiator the Gammacell 3000 and the BioBeam 8000 was appeared to 14% and 22%, respectively.

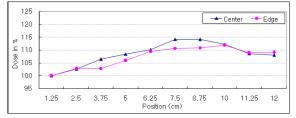


Fig. 1. The percent distributions in the Gammacell 3000 irradiator by using the glass dosimeter.


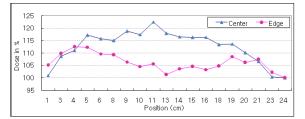


Fig. 2. The percent distributions in the BioBeam 8000 irradiator by using the glass dosimeter.

Fig. 3, Fig. 4 shows the results of 2-demenizional dose distribution having 2 mm spacing on the film image analyzed from spatial dose distribution using the Gafchromic EBT film for different doses.

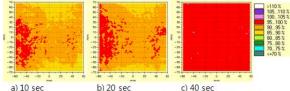

Dose variation ±7.5% Dose variation ±7.5% Dose variation ±2.5%

Fig. 3. The uniformity of dose distributions in the Gammacell 3000 irradiator by using the Gafchromic EBT film, a). 10 sec, b). 20 sec, and c). 40 sec.

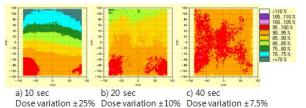


Fig. 4. The uniformity of dose distributions in the BioBeam 8000 irradiator by using the Gafchromic EBT film, a). 10 sec, b). 20 sec, and c). 40 sec.

### Conclusions

The measurement results of uniformity of delivered dose within the device showed variation more than 14% between middle point and the lowest position at central axis. Therefore, to maintain dose variation within 10%, the criteria of useful dose distribution, for research radiation effects, the irradiated specimen located at central axis of the container should be placed within 30 mm from top and bottom surface, respectively. In addition, for measurements using the film, the variations of dose distribution were more then 50% for the case of less than 10 second irradiation, mostly within 20% for the case of more than 20 second irradiation, respectively. Therefore, the irradiation experiments using the BioBeam 8000 irradiator are recommended to be used for specimen required at least more than 20 second irradiation time.

#### Reference

- Technical Report, Explanation Material of RPL Glass Dosimeter. Small Element System. Asahi Technev Glass Corporation, (2000)
- 2. Self-Developing film for radiotherapy dosimetry, GAFCHROMIC® EBT,(2007)