The Matrix Effect of Biological Concomitant Element on the Signal Intensity of Ge, As, and Se in Inductively Coupled Plasma/Mass Spectrometry

Kyung-Su Park,* Sun-Tae Kim, Young-Man Kim, Yunje Kim,* and Won Lee^{‡,*}

Advanced Analysis Center, Korea Institute of Science and Technology, Seoul 136-791, Korea Doping Control Center, Korea Institute of Science and Technology, Seoul 136-791, Korea Research Institute for Basic Sciences and Department of Chemistry, Kyunghee University, Seoul 130-701, Korea Received April 4, 2002

The non-spectroscopic interference effects that occurred in inductively coupled plasma/mass spectrometry were studied for Ge, As and Se in human urine and serum. Many biological samples contain Na, K, Cl and organic compounds, which may cause the enhancement and depression on the analyte signal. The effect of 1% concomitant elements such as N, Cl, S, P, C, Na, and K on a 100 µg/L germanium, arsenic and selenium signal has been investigated by ICP/MS. The interference effects were not in the same direction. It appeared that concomitant elements such as Cl, S, and C induce an enhancement effect, whereas N and P did not show any significant effect. And, Na and K caused a depression. We have found a link between the abundance of analytes and the ionization potential of concomitant elements (eV), except carbon and nitrogen.

Key Words: Matrix effect, Serum, Urine, ICP/MS

Introduction

Inductively coupled plasma (ICP) as an ionization source for mass spectrometry (MS) has gained general acceptance as a method to analyze the inorganic elements. The advantages of ICP/MS are well known; low limits of detection and simultaneous multi-element determination with the possibility of measuring isotopic abundance. The main limitations are the observation of isobaric interferences, mostly due to the formation of molecular compounds, in particular, oxides and non-spectroscopic interferences (matrix effect).^{1,2} The accuracy of As and Se determinations is well known to be seriously affected by the occurrence of spectral interferences.³⁻¹⁹ as the resolution of quadrupole mass spectrometers is not sufficient to resolve atomic and molecular species having the same nominal mass. But recently, spectral interferences has been overcome by collision cell or dynamic

Non-spectral interferences is usually defined as matrixinduced signal variations (both suppression and enhancement) and are therefore often referred to as matrix effects. F. Vanhaecke et al. have reported that matrix-induced signal variation is dependent on the analyte mass number.²⁰⁻²² We know from experience that the calibration of As and Se may be difficult, especially in biological and clinical matrices. These matrix effects come from the sample matrix (Na. K and organic compounds). However, there is a the number and variety of the mechanisms proposed to explain these effects. 23-26

In the present work, we report some preliminary investigations on the matrix effect by CH₃OH, C₂H₃OH and mannitol as the organic matrix and NaOH, KOH as the alkali

matrix. A study of non-spectral interferences in the determination of Ge. As and Se is described.

Experimental Section

Instrumentation. The ICP/MS instrument in this work was a Perkin-Elmer SCIEX Elan 5000 (Norwal, Connecticut, U.S.A.) ran at normal resolution. A cross-flow nebulizer associated with a double-pass spray chamber was used. The operating parameters were power of 1.0 kW, outer gas flow rate of 15 L/min., auxiliary 0.85 L/min, and carrier gas flow rate of 0.92 L/min. Mass flow controllers were used to measure flow rates. The mass spectrometer was set to sample ion intensities (peak jumping) at the analyte masses m/z 70, 72, 73 and 74 for Ge⁻, m/z 75 for As⁻ and m/z 77, 78 and 82 for Se⁻. A peristaltic pump was used for conventional sample introduction. Instrument sensitivity was optimized by varying on instrumental settings, using aqueous standard solutions of Pb. Rh and Mg. The instrument adjustments included the physical positioning of the plasma relative to the mass spectrometer, the ion lens voltages and the R. F. power relative to the argon plasma. Once optimized for germanium, arsenic and selenium, these variables were kept constant throughout all experiments involving the same element. The aerosol carrier gas flow was adjusted for each analyte until optimum signal intensity was reached. The ICP/MS setting is given in Table 1.

Materials and chemicals. Scrum and urine samples were obtained from 30 healthy individuals living in the Seoul region. Standard solutions were prepared from 10 mg/L stock solutions of Ge. As and Se (Spex Industries Inc., Edison, New Jersey, U.S.A.). All standard solutions were made up in 1% HNO₃. High purity reagents were used and the water was produced in a Millipore Super-Q apparatus (Millipore, Milford, MA, U.S.A.).

^{*}Corresponding Authors, Prof. Won Lee(Fax: +82-2-968-6584; e-mail: wonlee \(\hat{a}^2\) khu.ac.kr), Dr. Kvung Su Park(e-mail: pks/\(\hat{a}^2\) kist.re.kr)

Table 1. Instrumental and analytical parameters for ICP/MS

| Descriptions | Conditions |
|-----------------------------|-------------------------------------|
| R.F. generator | Free-running type, 40 MHz |
| R.F. power | 1000 W |
| Induction coil | 3-turn, 1 8 in. copper, 2.6 mm i.d. |
| Sampling depth | 7 mm from load coil, on center |
| Coolant gas flow rate | 15.0 L min. |
| Auxiliary gas flow rate | 0.85 L min. |
| Nebulizer gas flow rate | 0.92 L min. |
| Mass flow controller | 4-channel |
| Sample introduction | Peristaltic pump |
| Sample uptake flow | 1.0 mL/min. |
| Nebulizer | Cross-flow type |
| Spray chamber | Double pass type (Scott type) |
| Torch | Demountable |
| Interface cones | Nickel |
| Mass analyzer | Quadrupole |
| Vacuum system | Turbo molecular pumps |
| Quadrupole chamber | 5×10^{-7} torr |
| Quantitative mode | |
| Replicate time (ms) | 300 |
| Dwell time (ms) | 100 |
| Sweeps/Reading | 3 |
| Reading/Replicate | 1 |
| Number of Replicates | 5 |
| Points Spectral peak | 3 |
| Ge Mass | 70, 72, 73, 74 |
| As/Mass | 75 |
| Sc/Mass | 77. 78. 82 |
| Sean mode | Peak hopping |
| Resolution | 0.9 |
| Total Quant Mode Mass Range | 67-87 |

Sample preparation. 1 mL of a serum or urine samples were digested with 5 mL c-HNO₃. After digestion, the solutions were quantitatively transferred into calibrated flasks (10 mL) and made up with distilled, de-ionized water.

Five replicates of each sample were prepared.

Standard solutions and calibration solutions were prepared from the stock solutions as required.

Procedure. The signal intensities obtained from the standard solutions of Ge, As and Se were measured after conventional aspiration into the ICP/MS instrument. Before introduction into ICP/MS. Ge, As and Se were each diluted by a different matrix solution. The standard solution of 100 μ g/L Ge, As and Se in 1% matrices was used for the matrix effect on analyte signal. We investigated matrix effects on the signal for Ge, As and Se in the presence of increasing levels of matrices (0. 1000, 5000, 10000, 50000 μ g/L).

Results and Discussion

Selection of mass for analyte. Noticeable interference was found for Ge. As and Se, showing that below m/z = 85 the background spectrum is not free of interference peaks. The background spectrum is not free of interference peaks. The background spectrum is not free of interference peaks. The background spectrum is not free of interference peaks. The background spectrum is not free from oxides or argon dimer ion (40 Ar 40 Ar) has been reported. The background so the isotope m/z = 77 and 82 was chosen instead of m/z = 80 for Se. Because the m/z 75 (40 Ar 35 Cl) was a possible interfering ion, it was compensated for by the m/z 77 (40 Ar 37 Cl).

Sensitivity of individual metals by concomitant elements. Acid digested urine and serum samples contain Na, K. Cl and organic compounds (Table 2, 3). These interfering matrices can be an enhancer or suppress the signal intensities of the Ge, As and Se ions in ICP/MS. The ICP/MS signal intensities of three metals dissolved in different matrixes is shown in Figure 1. The signal intensity (sensitivity) of Ge, As and Se in ICP/MS was recorded after conventional aspiration of aqueous solutions of the standard substances. The results show that concomitant elements such as Cl, S, and C, induce an enhancement effect, whereas N and P did not have any significant effect. On the other hand, Na and K caused a suppression.

The signal intensity is enhanced by 77% (4 Ge). 210% (5 As) and 156% (8 ²Se) in 1% HCl matrix, and by 26%

Table 2. The concentrations of the concomitant elements in digested serum

| | Mf. o | | | | | |
|--------|---------|---------|----------|--------------------|--------------------|-------------------------|
| | Na | K | Р | S | С | CI |
| Mean | 0.0240 | 0.0221 | 0.0107 | ⊴0.1 | 0.14 | 0.057 |
| S.D. | 0.00086 | 0.00166 | 0.000833 | 0 | 0.0296 | 0.0233 |
| Method | AAS | AAS | ICP/AES | Elemental Analysis | Elemental Analysis | Potentiometric Titrator |

N=10

Table 3. The concentrations of the concomitant elements in digested urine

| | Wt.º o | | | | | |
|--------|--------|--------|---------|--------------------|--------------------|-------------------------|
| | Na | K | P | S | С | Cl |
| Mean | 0.0676 | 0.050 | 0.0073 | . 0.1 | 0.0410 | 0.0229 |
| S.D. | 0.0371 | 0.0936 | 0.00462 | 0 | 0.0619 | 0.0528 |
| Method | AAS | AAS | ICP/AES | Elemental Analysis | Elemental Analysis | Potentiometric Titrator |

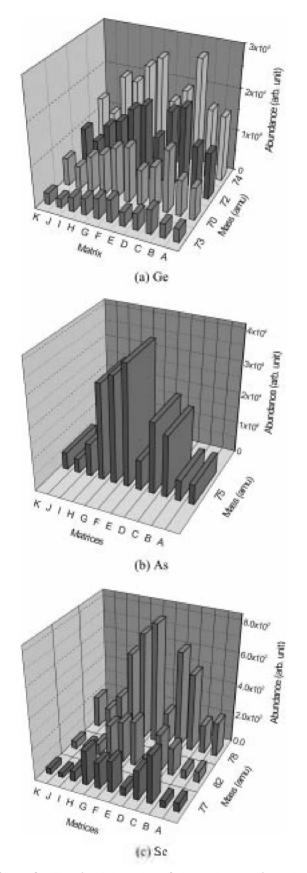


Figure 1. The signal intensity of three elements in various matrices: A: H₂O, B: IINO₃, C: IICl, D: H₂SO₄, E: II₃PO₄, F: CH₃OH, G: C₂H₅OH, H: Mannitol, I: NH₄OH, J: NaOH, K: KOH.

Table 4. The ratio of the intensity of ionic line and the intensity of atomic line for Pd in 1% matrices by ICP/AES^a

| | Ita | tensity of ionic l | ine | | |
|------------------|--------------------------|--------------------|----------|--|--|
| No. | Intensity of atomic line | | | | |
| _ | HNO ₃ | NaNO ₃ | Mannitol | | |
| Mean | 0.107 | 0.093 | 0.111 | | |
| S.D. | 0.000805 | 0.00326 | 0.000504 | | |
| Difference (° o) | _ | -12.9 | :3.36 | | |

^aConcentration of Pd is 10 μg/L, N=5

(14 Ge). 257% (75 As) and 251% (82 Se) in 1% H₂SO₄ matrix (Figure 1). In case of 1% HNO₃, the signals of these analytes showed no variation.

In the case of alkali metals (1% Na and K), the result shows the depressed signal intensity for analyte. Under these conditions the intensity of analytes was decreased by 8.8-6.7% (⁴Ge). 7.7-6.8% (⁵As) and 24.2-24.7% (⁸²Se) (Figure 1). And the emission intensity of ICP/AES for Pd (10 mg/L) standard solution shows a depression in 1% NaNO₃ of 12.9% and enhancement in 1% mannitol of 3.3% (Table 4). We projected the same effects of the matrix for Pd for Ge, As and Se in ICP/MS. From these results, we theorized that the lower sensitivity of analyte resulted from the incomplete ionization in the plasma and/or the recombination of Ge, As and Se with electrons.

The results show that the maximum metal signal enhancement is achieved when 1% organic matrices is added to the analyte solution. The signal intensities are enhanced by 24-67% (Ge). 383-501% (As) and 191-289% (Se) (Figure 1). The signal intensity of As dissolved in mannitol solution at the concentration of 1% as C increased to 383% compared with aqueous solution. This increase was similar to 1% methanol and ethanol.

The effects of matrix concentrations are summarized in Figure 2. The signal intensity of Ge, As and Se can increase 1.5%. 20% and 11%, respectively, for a CH₃OH concentration of only 0.1% as C. And the signal intensity in 0.1% $NaNO_3$ as Na decreased bt more than 10%.

Allain *et al.*, selected Eu as an internal standard.³¹ The use of Eu is not adequate with most elements because the interference effects are not in the same direction (Figure 4).

Variation of signal intensity of analyte by ionization potential of matrix. We have investigated the influence of ionization energy, because the enhancement effect would be linked to the chemical or physical properties of the elements. A plot of the enhancement percentage for matrices at 1% as a function of the ionization energy for each element is given in Figure 3. Because nitrogen in the anunonium or nitrate did not show the enhancement or depression effect, their salts were used. The results show that the signal intensity of analyte is not completely dependent on the ionization potential of matrix. But there would be an effect between signal intensity and ionization potential (eV) of concomitant element. Carbon (11.20 eV) is the only noticeable exception in this range. It should be noted that nitrogen (14.53 eV) also does not exhibit an enhancement effect for sensitivity.

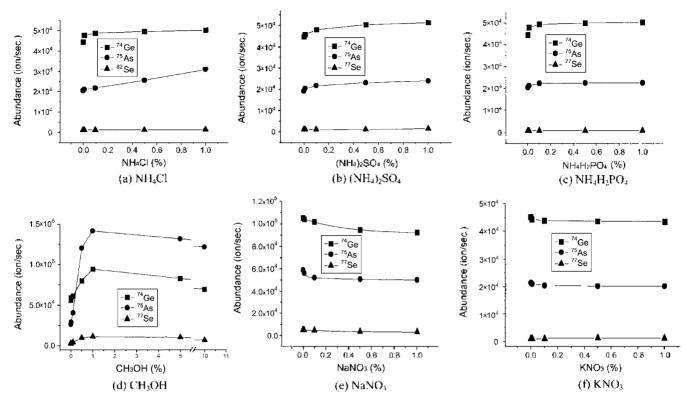


Figure 2. The effect of some matrix concentrations on the abundance of Ge, As and Se.

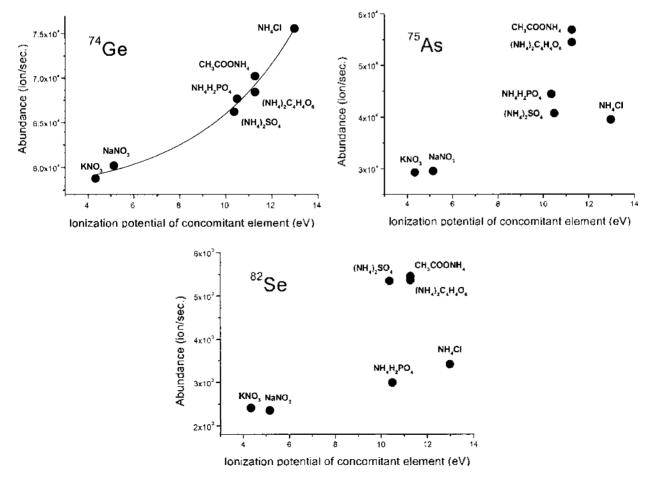


Figure 3. The effect of the concomitant elements on the abundance of Ge, As and Se.

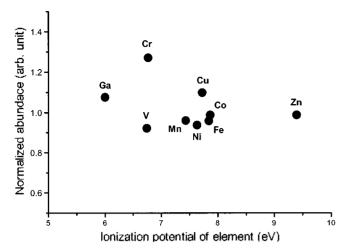


Figure 4. Normalized abundance of elements in 1% CH₃OH.

Conclusions

This work shows that interfering matrices such as Na. K and Carbon can cause signal enhancement or depression. Biological samples were digested in acid media containing concomitant elements in the range of 0.01%-0.18%. The enhancement of As can reach 20% for a CH₃OH concentration of only 0.1% as C. There is an effect between the abundance and ionization potential of the concomitant elements (eV), except carbon and nitrogen.

Acknowledgment. This work was supported by grant No. 1999-1-124-001-3 from the interdisciplinary research program of the KOSEF.

References

- 1. Tan. S. H.: Horlick, G. Appl. Spectrosc. 1988, 40, 445.
- Vaughan, M. A.; Horlick, G. Appl. Spectrosc. 1986, 40, 434.
- 3. Gray, A. L. Spectrochim, Acta, Part B 1986, 41, 151.
- Date, A. R.; Cheung, Y. Y.; Stuart, M. E. Spectrochim. Acta. Part B 1987, 42, 3.

- Beauchemin, D.; Craig, J. M. Spectrochim. Acta, Part B 1991, 46, 603
- McLaren, J. W.; Beauchemin, D.; Berman, S. S. J. Anal. At. Spectrom. 1987, 2, 277.
- Alves, L. C.: Wiederin, D. R.; Houk, R. S. Anal. Chem. 1992, 64, 1164
- Branch, S.; Ebdon, L.; Ford, M.; Foulkes, M.; Oneill, P. J. Anal. At. Spectrom. 1991, 6, 151.
- Branch, S.; Corns, W. R.; Ebdon, L.; Hill, S.; Oneill, P. J. Anal. At. Spectrom. 1991, 6, 155.
- Buckley, W. T.: Budae, J. J.: Godfrey, D. V.: Koenig, K. M. Anal. Chem. 1992, 64, 299.
- 11. Evans, E. H.; Ebdon, L. J. Anal. At. Spectrom. 1989, 4, 299.
- 12. Evans, E. H.; Ebdon, L. J. Anal. At. Spectrom. 1990, 5, 425.
- Hansen, S. H.; Larsen, E. H.; Pritzl, G.; Comett, C. J. Anal. At. Spectrom. 1992, 7, 629.
- Lyon, T. D. B.: Fell. G. S.: Hutton, R. C.: Faton, A. N. J. Anal. At. Spectrom. 1988, 3, 265.
- Lyon, T. D. B.; Fell, G. S.; Hutton, R. C.; Eaton, A. N. J. Anal. At. Spectron. 1988, 3, 601.
- Munro, S.; Ebdon, L.; McWeeny, D. J. J. Anal. At. Spectron. 1986, 1, 211.
- Ridout, P. S.: Jones, H. R.; Williams, J. G. Analyst 1988, 113, 1383.
- Sheppard, B. S.; Caruso, J. A.; Heitkemper, D. T.; Wolnik, K. A. *Analyst* 1992, 117, 971.
- 19. Tan. S. H.: Horlick, G. Appl. Spectrosc. 1986, 40, 445.
- Vanhaecke, F.: Vanhoe, H.: Dams, R.: Vandecasteele, C. Talanta 1992, 39, 373.
- Thompson, J. J.; Houk, R. S. Appl. Spectrosc. 1987, 41, 801.
- 22. Doherty, W. Spectrochim, Acta, Part B 1989, 44, 263.
- Beauchemin, D.; Mclaren, J. W.; Berman, S. S. Spectrochim, Acta. Part B 1987, 42, 467.
- 24. Olivares, J. A.: Houk, R. S. Anal, Chem. 1986, 58, 20.
- 25. Gregoire, D. C. Spectrochim, Acta, Part B 1987, 42, 895.
- Gilson, G. R.; Douglas, D. J.; Fulford, J. E.; Halligan, K. W.; Tanner, S. D. Anal. Chem. 1988, 60, 1472.
- 27. Evans, E. H.; Giglio, J. J. J. Anal. At. Spectrom. 1993, 8, 1.
- Vanhoe, H.; Goossens, J.; Moens, L.; Dams, R. J. Anal. At. Spectrom. 1994, 9, 177.
- Munro, S.; Ebdon, L.; McWeeny, D. J. J. Anal. At. Spectrom 1986, J. 211.
- Amarasiriwardena, C. J.; Lupoli, N.; Korrick, S.; Hu, H.; Potula, V. Analyst 1998, 123, 441.
- Abdallah, M. H.; Mermet, J. M. J. Quant. Spectrosc. Radiat. Transfer 1978, 19, 83.