



ARTICLE

## Multi-Residue Analysis of Fipronil and Its Metabolites in Eggs by SinChERS-Based UHPLC-MS/MS

Keguang Han<sup>1,†</sup>, Jin Hua<sup>2,†</sup>, Qi Zhang<sup>3,†</sup>, Yuanhui Gao<sup>2</sup>, Xiaolin Liu<sup>2</sup>, Jing Cao<sup>4</sup>, and Nairui Huo<sup>1,\*</sup>

<sup>1</sup>College of Veterinary Medicine, Shanxi Agricultural University, Taigu 030801, China

<sup>2</sup>Taiyuan Customs of the People's Republic of China, Taiyuan 030024, China

<sup>3</sup>College of Food Science, Southwest University, Chongqing 400715, China

<sup>4</sup>Beijing Vocational College of Agriculture, Beijing 102442, China

OPEN ACCESS

Received June 8, 2020

Revised August 23, 2020

Accepted August 25, 2020

\*Corresponding author : Nairui Huo  
College of Veterinary Medicine at Shanxi  
Agricultural University, Taigu 030801, China  
Tel: +86-13935452616  
Fax: +86-03546288160  
E-mail: [tnrhuo@163.com](mailto:tnrhuo@163.com)

\*ORCID  
Keguang Han  
<https://orcid.org/0000-0002-4596-6226>  
Jin Hua  
<https://orcid.org/0000-0001-9514-1487>  
Qi Zhang  
<https://orcid.org/0000-0002-9021-3960>  
Yuanhui Gao  
<https://orcid.org/0000-0001-9418-0329>  
Xiaolin Liu  
<https://orcid.org/0000-0002-4785-0657>  
Jing Cao  
<https://orcid.org/0000-0001-8352-4923>  
Nairui Huo  
<https://orcid.org/0000-0001-6046-0699>

† These authors equally contributed and should be regarded as co-first authors.

**Abstract** A method for simultaneous detection of fipronil (F) and its metabolites fipronil desulfinyl (FD), fipronil sulfide (FS), fipronil sulfone (FSO) in chicken eggs was applied and validated. It includes single-step, cheap, effective, rugged, safe-based method (SinChERS) for sample preparation and ultra high performance liquid chromatography coupled with mass spectrometry (UHPLC-MS/MS) for chemical analysis. Results suggested that formic acid enhanced the recovery of 4 target residues and 1% supplementation to acetonitrile gained higher recoveries than that of 5%. SinChERS integrated extraction and clean-up steps into one, with shorter time (1.5 h) to operate and higher recoveries (97%–100%) than HLB, Envi-Carb-NH<sub>2</sub> and quik-easy-cheap-effective-rugged-safe method (QuEChERS), and it consumed the smallest volume of extracting solvent (10 mL) as QuEChERS. Quantitative analyses using external standard method suggested the linear ranges of 4 target compounds were 1–20 µg/L with R<sup>2</sup>>0.9947. The limit of detection (S/N>3) and quantification (S/N>10) were 0.3 µg/kg and 1 µg/kg. Recoveries ranged from 89.0% to 104.4%, and the relative standard deviations (n=6) at 1, 10, and 20 µg/kg were lower than 6.03%. Thirty batches of domestic eggs (500 g each) were detected by the established SinChERS-based UHPLC-MS/MS and no target residues were detected in all samples. The method developed in this study is a rapid, sensitive, accurate and economic way for multi-residue analysis of fipronil and its metabolites in eggs.

**Keywords** single-step, cheap, effective, rugged, safe-based method (SinChERS), fipronil, egg, ultra high performance liquid chromatography coupled with mass spectrometry (UHPLC-MS/MS)

## Introduction

Eggs tainted with fipronil (F) were firstly reported in Belgium in 2017, then in the European Union, Korea, Hong Kong and Taiwan of China (Erickson, 2017). Soon afterwards chicken meat was also unveiled contaminated with fipronil (Stefanka et al.,

2017). Fipronil is a kind of broad-spectrum insecticide of phenylpyrazole group sprayed in poultry house to prevent and treat ectoparasite infestation (Cochran et al., 2015), leading to its bioaccumulation in eggs. Additionally, misuse or abuse of fipronil can also cause fipronil residue in eggs. Long-term low-dose or short-time high-dose intake of egg-born fipronil and its metabolites can put liver, thyroid, kidney and nervous system under health risks (Kitulagodage et al., 2011; Simon-Delso et al., 2015). The maximum residue limit (MRL) of fipronil and its metabolites in vegetables, fruits, grain, and oil have been stipulated in many related food standards and regulations. The MRL of fipronil and its metabolites in eggs stipulated by the Codex Alimentarius Commission (CAC) is 20 µg/kg, the European Union 5 µg/kg, the USA 30 µg/kg, while it hasn't been set in China.

Nowadays, various analytical methods have been introduced to determine fipronil and its metabolites in diverse matrices such as fruits (Duhan et al., 2015), vegetables (Kaur et al., 2015), peanut, soil et al. (Li et al., 2015), including gas chromatography (GC) (Guo et al., 2008), GC in tandem with mass spectrometry (GC-MS/MS) (Liu et al., 2019; Ramasubramanian et al., 2014), high performance liquid chromatography (HPLC) (Liu et al., 2008; Neagu et al., 2015), and liquid chromatography in tandem with mass spectrometry (LC-MS/MS) (Raju et al., 2016). Besides, QuEChERS-based method was also applied to determine fipronil in many matrices including eggs (Sack, et al., 2011; Xia et al., 2010), and it was popularized worldwide and proved effective, sensitive and accurate (Aruna et al., 2015; Yu et al., 2015).

The method of single-step, cheap, effective, rugged, safe-based method (SinChERS) is characterized by single-step, cheap, effective, rugged and safe. The validation parameters are based on the standards of AOAC Official Method 2007.01 and EN15662. It is a novel and proprietary way for sample preparation developed by the Anybond Technologies, Tianjin, China. QuEChERS was invented in 2003 by American chemists and widely used as a sample preparation technology in pesticides residue detection. Multi-wall carbon nanotubes (MWCNTs) was used in SinChERS to absorb different interfering co-extracts, while QuEChERS use primary secondary amine (PSA). Compared with PSA, MWCNTs have larger specific surface area and widely distributed mesopores and micropores in structure and thus exhibited excellent adsorption ability than PSA.

Animal-derived food matrix is so chemically complicated that sample preparation is extremely important for trace analysis, during which how to avoid the interference caused by co-extraction of non-target substances is the biggest challenge and hurdle. Sample preparation always accounts for 2/3 time of the whole process of analysis, directly impacting the efficiency and the accuracy of the quantification analysis and the instrument performance. The better method for sample treatment, the higher efficiency and precision of trace analysis.

Different from QuEChERS, SinChERS integrates extraction and purification steps into one single operation, avoiding loss of the target analyte during solvent transferring. In this one-step operation, only 5 mL extraction solvent was consumed for 10 g egg matrix and purification effect was good enough for later ultra high performance liquid chromatography coupled with mass spectrometry (UHPLC-MS/MS) analysis. It possess a higher recovery rate than QuEChERS and dispersive solid phase extraction (d-SPE).

In this study, SinChERS was applied for sample preparation and comparisons were made to QuEChERS and d-SPE in terms of time needed, solvent volume consumed, recovery rate and matrix effect (%ME) in the determination of fipronil and its metabolites fipronil desulfinyl (FD), fipronil sulfide (FS) and fipronil sulfone (FSO) by UHPLC-MS.

## Materials and Methods

### Chemicals and instrumentations

Standard reagents of F, FD, FS, and FSO (HPLC-grade, ≥99%) were supplied by Dr. Ehrenstorfer (Augsburg, Germany).

HPLC-grade solvents ammonium acetate, methyl alcohol, formic acid and acetonitrile were procured from Fisher Scientific (Fair Lawn, NJ, USA). Anhydrous magnesium sulfate ( $\text{MgSO}_4$ ) of analytical-grade was the product of Agela Technologies (Tianjin, China). Water was collected from a Milli Q purification system (Millipore, Molsheim, France).

Waters Oasis HLB Solid phase extraction (SPE) purification column (6 mL/500 mg) was provided by the Troody Analytical Instrument (Shanghai, China). Supelco Envi-Carb/LC- $\text{NH}_2$  purification column (6 mL/500 mg) was purchased from the Kanglin Instrument (Beijing, China). QuEChERS purification column was from Sepax Technologies (Guangzhou, China). SinChERS purification column was the product of Anybond Technologies (Tianjin, China, <http://www.anybond.com.cn/SinCHERS2>), and was filled with 900 mg  $\text{Na}_2\text{SO}_4$ , 150 mg MWCNTs, and 150 mg C18.

The UHPLC-MS/MS system consists of a Waters UPLC system (USA), a Quattro Premier XE quadrupole mass spectrometer (The Science of Waters, USA), and a electrospray ionization (ESI) interface source. Other instruments include ultrasonic instrument (KQ5200DE, Kun Shan Ultrasonic Instruments, China), advanced vortex mixer (EFO945601, TALBOYS, USA) and bath-typed nitrogen evaporator (OA-SYS, Organomation Associates, USA).

### Preparation of standard solutions

Stock solutions of standards F, FD, FS, and FSO were prepared by dissolving each into acetonitrile to 1,000 mg/L, stored in dark vials at  $-18^\circ\text{C}$ . Working standard solution was a mixture of 4 standards, 1  $\mu\text{g}/\text{mL}$  each, diluted with acetonitrile. The calibration standards (1, 2, 5, 10, and 20  $\mu\text{g}/\text{L}$ ) were prepared by stepwise dilution of the above working solution of each analyte (1  $\mu\text{g}/\text{mL}$ ) with blank matrix solution.

### Sample preparation

The egg content was collected and homogenized as egg sample. Five gram of egg sample along with 10 mL 1% formic acid acetonitrile was added to a 50 mL centrifuge tube and blended for 3 minutes followed by 10 minute ultrasonic treatment. Then 2 g anhydrous  $\text{MgSO}_4$  was added and mixed thoroughly. After centrifuging at 4,000 r/min for 5 minutes, a SinChERS purification column was vertically placed into the centrifuge tube and manually pressed downwards to the tube bottom. With the movement of the purification column, organic phase of the sample entered into the reservoir tank of the column through the purification bed. Approximately 2 mL liquid in the reservoir tank was transferred into a new centrifuge tube and dried at  $40^\circ\text{C}$  in nitrogen and then dissolved into 1 mL mobile phase containing 65% 1 mM ammonium acetate, and finally filtered through a 0.22  $\mu\text{m}$  syringe filter. The filtrate collected was the final sample ready for UHPLC-MS/MS analysis.

### Ultra high performance liquid chromatography (UHPLC) conditions

The injection volume was 10  $\mu\text{L}$ . A Shim-pack GIST C18 column ( $50 \times 2.1$  mm, 2  $\mu\text{m}$ , Shimadzu) was applied and the temperature was controlled at  $40^\circ\text{C}$ . The mobile phase was the mixture of 1 mM ammonium acetate (A) and methyl alcohol (B). The flow rate was set to 0.3 mL/min. Gradient elution programs were given in Table 1.

### Mass spectrometry (MS/MS) conditions

MS/MS analyses were conducted on a LCMS-8050 equipped with an ESI source. All analytes were scanned by triple quadrupole multiple reaction monitoring mode (MRM) of negative ionization. Temperatures of desolvation line (DL), heat block and interface were  $250^\circ\text{C}$ ,  $400^\circ\text{C}$  and  $300^\circ\text{C}$ , individually. The flow rate of nebulizer gas ( $\text{N}_2$ ), heater gas ( $\text{N}_2$ ) and

**Table 1. Combinations of the mobile phase and the corresponding elution time for UHPLC**

A (%)	B (%)	Time (min)
65	35	0
45	55	1.5
15	85	3.5
65	35	3.5
65	35	4.0

A (1 mM ammonium acetate) and B (methyl alcohol) are solutions making up the mobile phase. UHPLC, ultra high performance liquid chromatography.

drying gas (N<sub>2</sub>) was 3 L/min, 10 L/min and 10 L/min, respectively. Argon (Ar) was used as the collision gas. Detailed parameters for MRM transitions were listed in Table 2 and MRM chromatograms of 4 standards were showed in Fig. 1.

### Evaluation of matrix effect

%ME was defined by the following equation where A is the peak area of the standard solution dissolved in acetonitrile, B is the peak area of the standard solution dissolved in blank egg matrix (Choi et al., 2015). If ME% falls into the scope of -20 to 20, it indicates no ME; ME% < -20 indicates ion suppression; ME% > 20 indicates a signal enhancement.

$$\text{ME\%} = (B - A) / A \times 100$$

## Results and Discussion


### Selection of extraction solvent

Fipronil is a weak-polar compound that can be easily dissolved into organic reagents (Jacob et al., 2015). Acetonitrile is a routine extracting solvent in precipitating protein and other conjugates (Li et al., 2016). Formic acid, together with acetonitrile, can improve the recovery rate of polar compounds and other impurities (Liu et al., 2014).

Organic-grade egg liquid (negative samples) spiked with 5.0 µg/L working standard solution was used to select the extraction solvent for SinChERS. Acetonitrile, 1% formic acid acetonitrile (acetonitrile+1% formic acid) and 5% formic acid


**Table 2. MRM parameters in MS/MS for fipronil and its metabolites**

Pesticide	Molecular formula	Retention time (min)	Precursor ions (m/z)	Product ions (m/z)	Cone voltage (V)	Collision energy (eV)
Fipronil (F)	C <sub>12</sub> H <sub>4</sub> C <sub>12</sub> F <sub>6</sub> N <sub>4</sub> OS	1.96	436.50	330.50	20	15
				250.20	20	30
Fipronil desulfinyl (FD)	C <sub>12</sub> H <sub>4</sub> C <sub>12</sub> F <sub>6</sub> N <sub>4</sub>	2.76	387.29	351.29	25	17
				282.31	25	26
Fipronil sulfide (FS)	C <sub>12</sub> H <sub>4</sub> C <sub>12</sub> F <sub>6</sub> N <sub>4</sub> S	3.69	418.91	262.35	20	15
				383.57	20	28
Fipronil sulfone (FSO)	C <sub>12</sub> H <sub>4</sub> Cl <sub>2</sub> F <sub>6</sub> N <sub>4</sub> O <sub>2</sub> S	2.85	450.29	415.31	22	27
				282.27	22	17

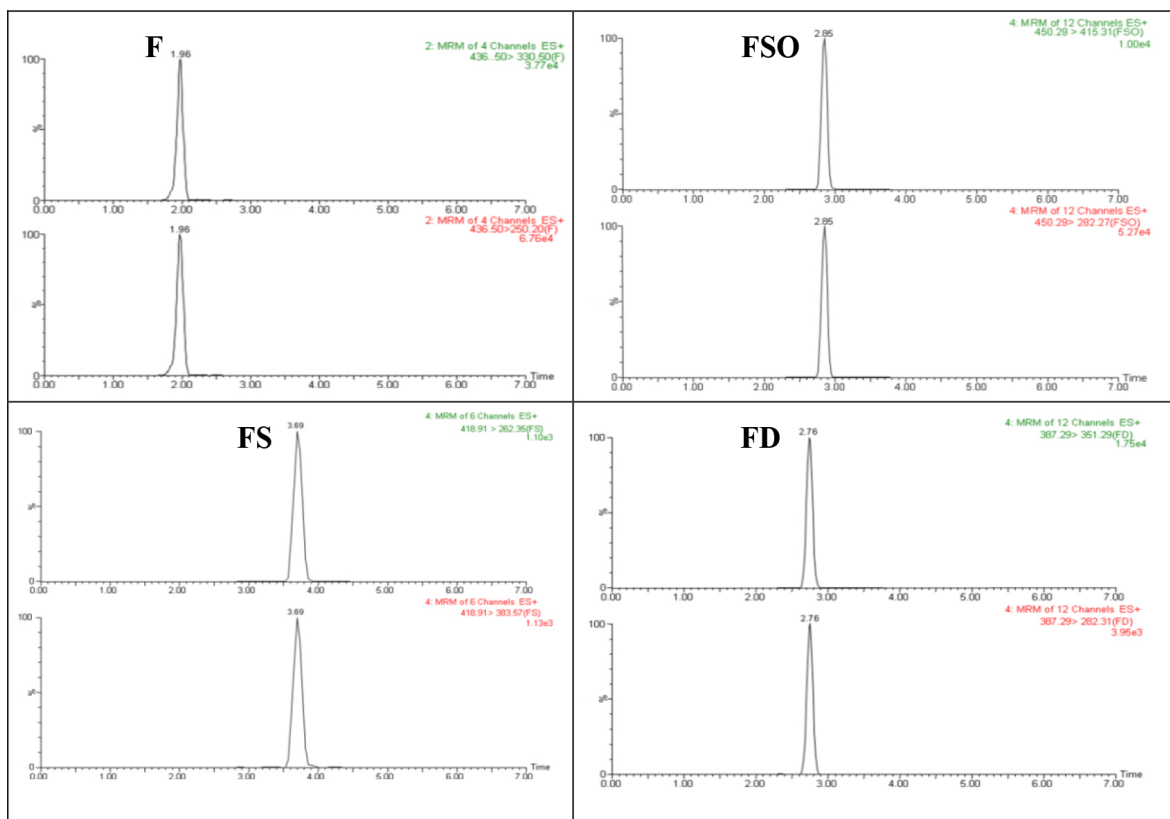

**Fig. 1.** Multiple reaction monitoring (MRM) chromatograms of F, FD, FS, and FSO standards. F: fipronil; FD, fipronil desulfanyl; FS, fipronil sulfide; FSO, fipronil sulfone.

acetonitrile (acetonitrile+5% formic acid) were recommended by the SinChERS manufacture as extraction solvents. Recoveries for the 4 analytes extracted by the above three tested candidate solvents were compared. It can be seen from Fig. 2 that all the recoveries were acceptable, ranging from 60% to 100%, and 1% formic acid acetonitrile showed the highest recovery rate for 4 target analytes.

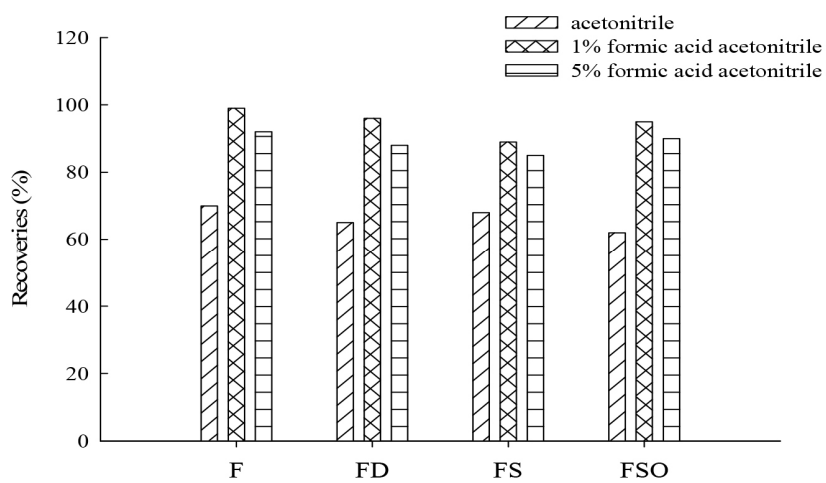

**Fig. 2.** Extract effects of different solvents for F, FD, FS, and FSO. F: fipronil; FD, fipronil desulfanyl; FS, fipronil sulfide; FSO, fipronil sulfone.

## Optimization of cleanup

Egg liquid contains a lot of fat, protein and some fat-soluble impurities (Jain et al., 2017; Stoddard et al., 2017). Different cleanup methods for fipronil and its metabolites in egg matrix have been reported, including QuEChERS (Shi et al., 2017) and d-SPE (Guo et al., 2017; Zhang et al., 2016). The prescribed d-SPE columns in Chinese trade standard SN/T 4039-2014 and the Chinese national standard GB 23200.34-2016 were HLB and Envi-Carb-NH<sub>2</sub>. Therefore, the cleanup effect of SinChERS was compared to QuEChERS, HLB, and Envi-Carb-NH<sub>2</sub> in this study in terms of recoveries, consumption of organic solvents and time needed, using matrix-matching external standard method in which 5.0 µg/L mixture of working standard solution was spiked into egg samples.

As shown in Table 3, different columns received acceptable recoveries from 69% to 100%. Generally, SinChERS>QuEChERS>HLB>Envi-Carb-NH<sub>2</sub>. SinChERS and QuEChERS consumed the least organic solvent (10 mL) and needed less time than two SPE columns, HLB and Envi-Carb-NH<sub>2</sub>. It's obvious that SinChERS clean-up column was even more time-saving than QuEChERS. Therefore, by comparison, SinChERS was a solvent-economic and time-saving option to extract the target analytes out of the complicated, fatty egg matrix, obtaining the best purification effect and the highest recoveries.

In similar studies, although the methods proved to be sensitive and reliable, the sample treatment is complex. In these studies, extraction and purification were separated operations, before purification, fipronil and its metabolites need to be salted out at -20°C (Guo et al., 2018) or evaporated to dryness (Liu et al., 2019) and be dissolved again. SinChERS avoids these steps and the operation is simplified, only ultrasonication and centrifugation are needed, thus make it convenient and time-saving. This novel rapid single-step extraction and cleanup method was applied to analyze up to 47 representative pesticide residues in vegetable and sauce products coupled with LC-MS/MS and GC-MS/MS detection, with modified cartridge fitted with multi-walled carbon nanotubes (MWCNTs), along with PSA and salts (Song et al., 2019).

## Matrix effect

Ionization suppression/enhancement could be brought about by sample matrix, sample preparation procedure or ionization types, leading to enhanced or weakened analytical signals and thus affecting the sensitivity and the precision of quantification.

**Table 3.** Comparison of 4 types of clean-up columns in terms of recovery, matrix effect, solvents and time consumption

Cleanup methods		HLB	Envi-Carb-NH <sub>2</sub>	QuEChERS	SinChERS
Recoveries (%)	F	76	70	92	100
	FD	85	73	94	99
	FS	77	68	89	88
	FSO	74	69	95	97
Matrix effect (%)	F	-1.00	-4.40	-2.70	1.10
	FD	17.60	7.60	11.40	10.60
	FS	20.40	6.90	16.30	16.20
	FSO	19.30	8.70	14.50	13.40
Solvent consumption		Acetonitrile 130 mL	Acetone 40 mL, dichloromethane 30 mL	Acetonitrile 10 mL	Acetonitrile 10 mL
Time spent (h)		4.5	6.5	2.0	1.5

F: fipronil; FD, fipronil desulfinyl; FS, fipronil sulfide; FSO, fipronil sulfone; QuEChERS, quik-easy-cheap-effective-rugged-safe method; SinChERS, single-step, cheap, effective, rugged, safe-based method.


ESI is more prone to incur such effects than atmospheric pressure chemical ionization, especially when other compounds are eluted together with the analyte of interest. Sample matrix is another principal factor especially in LC-MS/MS and usually exhibited ionization inhibiting effects on ESI.

Mean values of ME% were presented in Table 3. Except HLB column for FS extraction, all ME% values fell into the scope of -20 to 20, indicating egg matrix showed no effects on UHPLC-MS/MS signals when treated by these cleanup columns. By comparison, the ME of HLB column was the biggest followed by QuEChERS and SinChERS. Envi-Carb-NH<sub>2</sub> achieved the smallest values since it consisted of SPE adsorbent of Supelclean Envi-Carb (superstratum) and LC-NH<sub>2</sub> (substratum) (Wu et al., 2012). However, these two materials were quite more expensive than that of SinChERS or QuEChERS column. In overall, SinChERS was a suitable column for sample preparation for UHPLC-MS/MS in analyzing fipronil and its metabolites in complicated egg fluid.

## Method Validation

### Linearity

The linearity of the selected SinChERS method was evaluated using matrix-matched calibrations by spiking mixed standard solutions into blank egg samples to a final concentration of 1, 2, 5, 10, and 20 µg/L. Coefficient of determination ( $R^2$ ) and linear equations were obtained from the calibration curves (Fig. 3) drawn by plotting the peak areas against the

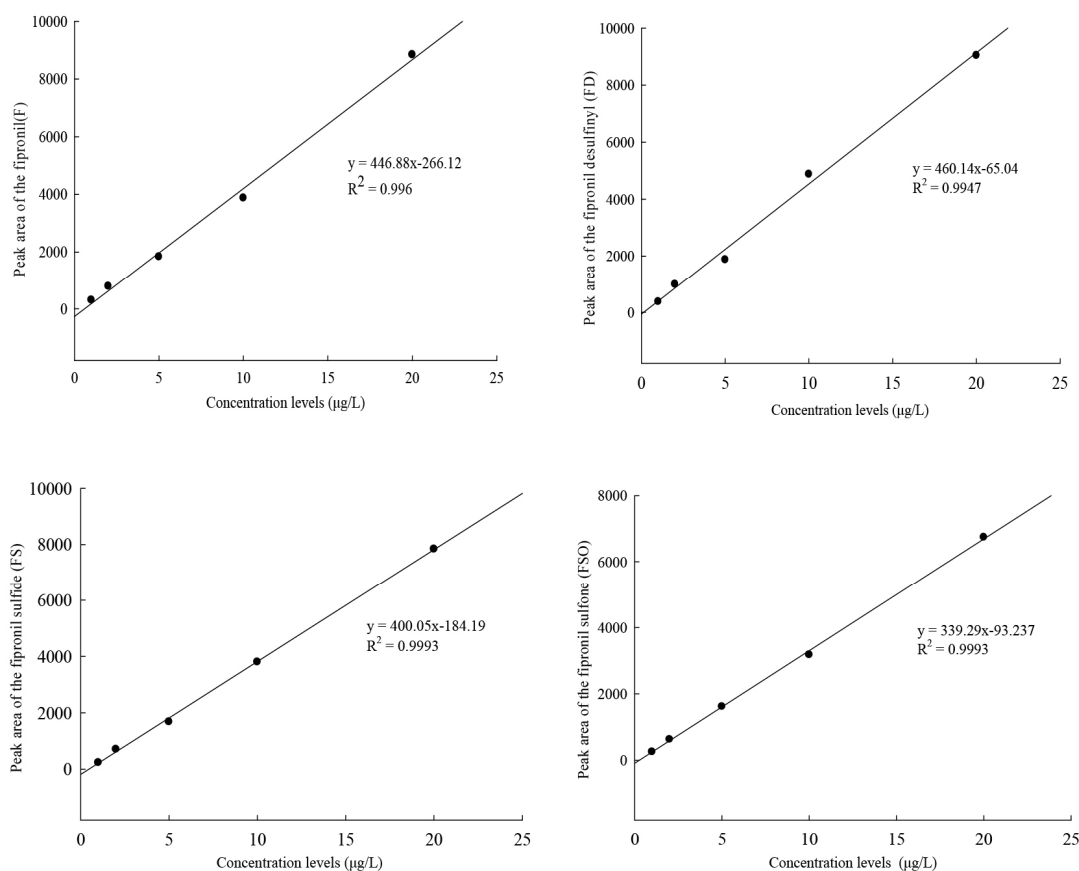

**Fig. 3.** Calibration curves of F, FD, FS, and FSO at 5 concentration levels. F: fipronil; FD, fipronil desulfinyl; FS, fipronil sulfide; FSO, fipronil sulfone.

concentrations of F, FD, FS, and FSO. As shown in Fig. 3, the calibration curves of fipronil and its metabolites exhibited sufficient linearity with  $R^2 > 0.9947$ , indicating that the proposed method could be applied for effective determination of trace pesticides in egg samples.

### Limit of detection and quantitation

To obtain the analytical limits of this method, the limit of detection (LOD) and the limit of quantitation (LOQ) were determined using signal-to-noise ratios (S/N) at 3 and 10, respectively. In the experiments described here, LODs and LOQs of all target analytes were 0.3  $\mu\text{g}/\text{kg}$  and 1  $\mu\text{g}/\text{kg}$ , meeting the criterion of Chinese GB 27417-2017.

### Recovery and precision

These were estimated by spiking the mixed standard solution to blank egg samples to the final level of 1  $\mu\text{g}/\text{kg}$ , 10  $\mu\text{g}/\text{kg}$ , and 20  $\mu\text{g}/\text{kg}$  in 6 duplicates according to the Chinese GB 2763-2016, which stipulated the MRL in most foods as 20  $\mu\text{g}/\text{kg}$ . Repeatability (intra-day precision) was assessed by recovery and precision was expressed as the relative standard deviations (RSDs).

As shown in Table 4, recoveries of F, FD, FS, and FSO were ranged from 89.00% to 104.40%. RSDs were 3.88%–6.03% at 1  $\mu\text{g}/\text{kg}$ , 4.15%–5.54% at 10  $\mu\text{g}/\text{kg}$ , and 2.22%–4.53% at 20  $\mu\text{g}/\text{kg}$ . All these recovery and RSD values met the criteria of

**Table 4.** Recoveries (%) and RSDs (%) of F, FD, FS, and FSO at 3 spiking levels by UHPLC-MS/MS (n=6)

Pesticide name	Spiked levels		1 $\mu\text{g}/\text{kg}$		10 $\mu\text{g}/\text{kg}$		20 $\mu\text{g}/\text{kg}$	
	Replicates	Measured value ( $\mu\text{g}/\text{kg}$ )	Average recoveries (%)	Measured value ( $\mu\text{g}/\text{kg}$ )	Average recoveries (%)	Measured value ( $\mu\text{g}/\text{kg}$ )	Average recoveries (%)	
Fipronil (F)	1	0.97	97.00	10.32	103.20	19.36	96.80	
	2	0.92	92.00	9.91	99.10	18.09	90.45	
	3	0.91	91.00	9.55	95.50	19.49	97.45	
	4	0.99	99.00	9.66	96.60	18.86	94.30	
	5	0.90	90.00	9.82	98.20	18.64	93.20	
	6	0.96	96.00	9.05	90.50	19.10	95.50	
	Mean value	0.94	94.17	9.72	97.18	18.92	94.62	
RSDs (n=6) (%)		3.88		4.34		2.72		
Recovery range (%)		90.00–99.00		90.50–103.20		90.45–97.45		
Fipronil desulfinyl (FD)	1	1.03	103.00	10.44	104.40	18.90	94.50	
	2	1.02	102.00	10.16	101.60	18.32	91.60	
	3	0.95	95.00	9.58	95.80	18.80	94.00	
	4	0.90	90.00	10.10	101.00	19.60	98.00	
	5	1.02	102.00	10.10	101.00	19.06	95.30	
	6	1.03	103.00	9.29	92.90	19.14	95.70	
	Mean value	0.99	99.17	9.95	99.45	18.97	94.85	
RSDs (n=6) (%)		5.46		4.27		2.22		
Recoveries range (%)		90.00–103.00		95.80–104.40		91.60–98.00		



**Table 4. Recoveries (%) and RSDs (%) of F, FD, FS, and FSO at 3 spiking levels by UHPLC-MS/MS (n=6) (continued)**

Pesticide name	Spiked levels		1 µg/kg		10 µg/kg		20 µg/kg	
	Replicates		Measured value (µg/kg)	Average recoveries (%)	Measured value (µg/kg)	Average recoveries (%)	Measured value (µg/kg)	Average recoveries (%)
Fipronil sulfide (FS)	1		0.91	91.00	9.39	93.90	20.00	100.00
	2		1.00	100.00	10.22	102.20	17.91	89.55
	3		0.90	90.00	10.01	100.10	18.58	92.90
	4		1.00	100.00	8.91	89.10	19.17	95.85
	5		0.91	91.00	9.64	96.40	18.83	94.15
	6		1.03	103.00	9.01	90.10	19.03	95.15
	Mean value		0.96	95.83	9.53	95.30	18.92	94.60
	RSDs (n=6) (%)		6.03		5.54		3.65	
Recoveries range (%)		90.00–103.00		89.10–102.20		89.55–100.00		
Fipronil sulfone (FSO)	1		0.94	94.00	10.25	102.50	17.94	89.70
	2		0.90	90.00	10.43	104.30	20.05	100.25
	3		0.98	98.00	9.58	95.80	19.47	97.35
	4		0.90	90.00	10.01	100.10	19.38	96.90
	5		0.90	90.00	9.58	95.80	19.62	98.10
	6		0.89	89.00	9.43	94.30	18.10	90.50
	Mean value		0.92	91.83	9.88	98.80	19.09	95.47
	RSDs (n=6) (%)		3.80		4.15		4.53	
Recoveries range (%)		89.00–98.00		94.30–104.30		89.70–100.25		

UHPLC-MS/MS, ultra high performance liquid chromatography coupled with mass spectrometry.

precision and accuracy. Thus the SinChERS, using 1% formic acid acetonitrile as extracting solvent, coupled with UHPLC-MS/MS was sensitive enough to detect and quantify part-per-billion level (ng/g or µg/kg) of pesticide residue in eggs, and it was suitable for multi-residue analysis of fipronil and its metabolites in eggs.

### Method application

The SinChERS-based UHPLC-MS/MS method was applied to real sample surveys. Thirty batches of domestic eggs were re-analyzed. These eggs have earlier been quarantined for other residues and all values were under the corresponding MRLs. They were collected from different farms of Shanxi province of China at different times. The results showed that neither fipronil nor its metabolites were detected above their LOQs.

### Conclusion

The SinChERS-based UHPLC-MS/MS method established in this study was sensitive for simultaneous analysis of fipronil and its metabolites in complicated egg matrix with high precision and reliability. The SinChERS, integrating sample extraction and clean-up steps into one, matching with 1% formic acid acetonitrile as extracting solvent, was faster, easier, more convenient, more solvent-economic and time-saving than HLB, Envi-Carb-NH<sub>2</sub> and QuEChERS columns. This sample

preparation procedure exhibited slight ME in later UHPLC-MS/MS analysis. In summary, the established method could play important roles in guaranteeing the safety of egg and egg products. Further studies needed to carry out to try the possibility of this method to analyze more other harmful residues simultaneously.

## Conflicts of interest

The authors declare no potential conflict of interest.

## Acknowledgments

This work was financially supported by the key R&D Programme of Shanxi Province (No. 201803D221023-4 and No. 201803D31223).

## Author Contributions

Conceptualization: Han K, Huo N. Data curation: Hua J, Zhang Q. Formal analysis: Gao Y. Methodology: Huo N, Hua J. Software: Liu X. Validation: Zhang Q, Cao J. Investigation: Han K, Hua J, Zhang Q. Writing - original draft: Han K, Gao Y. Writing - review & editing: Han K, Hua J, Zhang Q, Gao Y, Liu X, Cao J, Huo N.

## Ethics Approval

This article does not require IRB/IACUC approval because there are no human and animal participants.

## References

- Aruna M, Raju RM, Reddy KN. 2015. Method validation and residues of fipronil and its metabolites in citrus fruit and soil on gas chromatography- ECD by QuEChERS method. *J Chem Biol Phys Sci* 6:184-193.
- Choi S, Kim S, Shin JY, Kim MK, Kim JH. 2015. Development and verification for analysis of pesticides in eggs and egg products using QuEChERS and LC-MS/MS. *Food Chem* 173:1236-1242.
- Cochran RC, Yu L, Krieger RI, Ross JH. 2015. Postapplication fipronil exposure following use on pets. *J Toxicol Environ Health A* 78:1217-1226.
- Duhan A, Kumari B, Duhan S. 2015. Determination of residues of fipronil and its metabolites in cauliflower by using gas chromatography-tandem mass spectrometry. *Bull Environ Contam Toxicol* 94:260-266.
- Erickson BE. 2017. Illegal use of fipronil contaminates EU eggs. *C & En Global Enterprise* 95:21-21.
- Guo DH, Shi YY, Li Y, Yi XH, Deng XJ, Xiao WQ, Wang J, Li X, Liu H, Shen WJ. 2017. Rapid screening of fipronil and its metabolites in egg and egg products by solid phase extraction-liquid chromatography-quadrupole time-of-flight mass spectrometry. *J Chromatogr* 35:1216-1223.
- Guo M, Zong LG, Shan ZJ, Shi LL. 2008. Determination of fipronil in soils by accelerated solvent extraction and gas chromatography. *J Agro-Environ Sci* 27:2114-2117.
- Guo Q, Zhao S, Zhang J, Qi K, Du Z, Shao B. 2018. Determination of fipronil and its metabolites in chicken egg, muscle and

- cake by a modified QuEChERS method coupled with LC-MS/MS. *Food Addit Contam A* 35:1543-1552.
- Jacob CR, Soares HM, Nocelli RC, Malaspina O. 2015. Impact of fipronil on the mushroom bodies of the stingless bee *Scaptotrigona postica*. *Pest Manag Sci* 71:114-122.
- Jain S, Dhakal D, Anal AK. 2017. Bioprocessing of chicken meat and egg processing industries' waste to value-added proteins and peptides. In *Food processing by-products and their utilization*. Anal AK (ed). John Wiley & Sons, Hoboken, NJ, USA. pp 368-393.
- Kaur R, Mandal K, Kumar R, Singh B. 2015. Analytical method for determination of fipronil and its metabolites in vegetables using the QuEChERS method and gas chromatography/ mass spectrometry. *J AOAC Int* 98:464-471.
- Kitulagodage M, Buttemer WA, Astheimer LB. 2011. Adverse effects of fipronil on avian reproduction and development: Maternal transfer of fipronil to eggs in zebra finch *Taeniopygia guttata*, and *in ovo*, exposure in chickens *Gallus domesticus*. *Ecotoxicology* 20:653-660.
- Li M, Li P, Wang L, Feng M, Han L. 2015. Determination and dissipation of fipronil and its metabolites in peanut and soil. *J Agric Food Chem* 63:4435-4443.
- Li Y, Chen Z, Zhang R, Luo P, Zhou Y, Wen S, Ma M. 2016. Simultaneous determination of 42 pesticides and herbicides in chicken eggs by UHPLC-MS/MS and GC-MS using a QuEChERS-based procedure. *Chromatographia* 79:1165-1175.
- Liu D, Wang P, Zhu W, Gu X, Zhou W, Zhou Z. 2008. Enantioselective degradation of fipronil in Chinese cabbage (*Brassica pekinensis*). *Food Chem* 110:399-405.
- Liu H, Gong YX, Wang SY, Yu KY, Li JJ, Shen WJ, Ding T, Zhao ZY, Wu B, Zhang R, Yi XH, Deng XJ, Xu DM, Han S. 2019. Determination of fipronil and its metabolite residues in eggs and egg products by GC-MS/MS. *J Chin Mass Spectr Soc* 40:74-82.
- Liu W, Cheng S, Guo J. 2014. Anode modification with formic acid: A simple and effective method to improve the power generation of microbial fuel cells. *Appl Surf Sci* 320:281-286.
- Neagu M, Socenau G, Bucur A, Tanase A. 2015. Development and validation of a high performance liquid chromatography (HPLC) method for determination of fipronil and pyriproxyfen in a veterinary pharmaceutical product. *Int J Pharm Sci Drug Res* 7:123-125.
- Raju KS, Taneja I, Rashid M, Sonkar AS, Wahajuddin M, Singh SP. 2016. DBS-platform for biomonitoring and toxicokinetics of toxicants: Proof of concept using LC-MS/MS analysis of fipronil and its metabolites in blood. *Sci Rep* 6:22447.
- Ramasubramanian T, Paramasivam M, Jayanthi R, Chandrasekaran S. 2014. A simple and sensitive single-step method for gas chromatography-mass spectrometric determination of fipronil and its metabolites in sugarcane juice, jaggery and sugar. *Food Chem* 150:408-413.
- Sack C, Smoker M, Chamkasem N, Thompson R, Satterfield R, Masse C, Mercer G, Neuhaus B, Cassias I, Chang E, Lin Y, MacMahon S, Wong J, Zhang K, Smith RE. 2011. Collaborative validation of the QuEChERS procedure for the determination of pesticides in food by LC-MS/MS. *J Agric Food Chem* 59:6383-6411.
- Shi YY, Li Y, Yi XH, Deng X, Xiao W, Wang J, Li X, Liu H, Shen W. 2017. Determination of fipronil and its metabolites in egg and egg products by liquid chromatography-tandem mass spectrometry. *J Inspect Quar* 5:1-7.
- Simon-Delso N, Amaral-Rogers V, Belzunces LP, Bonmatin JM, Chagnon M, Downs C, Furlan L, Gibbons DW, Giorio C, Girolami V, Goulson D, Kreutzweiser DP, Krupke CH, Liess M, Long E, McField M, Mineau P, Mitchell EAD, Morrissey CA, Noome DA, Pisa L, Settele J, Stark JD, Tapparo A, Van Dyck H, Van Praagh J, Van der Sluijs JP, Whitehorn PR,

- Wiemers M. 2015. Systemic insecticides (neonicotinoids and fipronil): Trends, uses, mode of action and metabolites. *Environ Sci Pollut Res* 22:5-34.
- Song L, Han Y, Yang J, Qin Y, Zeng W, Xu S, Pan C. 2019. Rapid single-step cleanup method for analyzing 47 pesticide residues in pepper, chili peppers and its sauce product by high performance liquid and gas chromatography-tandem mass spectrometry. *Food Chem* 279:237-245.
- Stefanka B, Lubomir K, Gerhard B, Joerg S, Håkan E, John S, Piotr R, Hendrik E. 2017. Determination of the fipronil content in eggs. Joint Research Centre (JRC) Technical Reports. Publications Office of the European Union, Luxembourg.
- Stoddard MC, Yong EH, Akkaynak D, Sheard C, Tobias JA, Mahadevan L. 2017. Avian egg shape: Form, function, and evolution. *Science* 356:1249-1254.
- Wu N, Zhang YD, Liu W, Yi P, Wang ZF, Huang DS. 2012. The determination of eight pesticides in *Panax notoginseng* of Chinese traditional medicine using solid-phase extraction and gas chromatography-mass spectrometry. *Adv Mat Res* 485:68-71.
- Xia K, Atkins J, Foster C, Armbrust K. 2010. Analysis of cyromazine in poultry feed using the QuEChERS method coupled with LC-MS/MS. *J Agric Food Chem* 58:5945-5949.
- Yu Y, Xu D, Lu M, Zhou S, Peng T, Yue Z, Zhou Y. 2015. QuEChERS combined with online interference trapping LC-MS/MS method for the simultaneous determination of 20 polyfluoroalkane substances in dietary milk. *J Agric Food Chem* 63:4087-4095.
- Zhang M, Bian K, Zhou T, Song X, Liu Q, Meng C, He L. 2016. Determination of residual fipronil in chicken egg and muscle by LC-MS/MS. *J Chromatogr B* 1014:31-36.